

HOUSE LIGHT / STIMULUS LIGHT

ENV-215 / ENV-221

WIRING GUIDE

DOC-066

Rev. 1.2

Copyright ©2017

All Rights Reserved

Med Associates Inc.

P.O. Box 319

St. Albans, Vermont 05478

Phone: 802.527.2343

Fax: 802.527.5095


www.med-associates.com

CHAPTER 1 | WIRING

OPERATE CABLE

Connect the Operate cable any available Output on a Standard MED Connection Panel. By turning the corresponding output bit “ON” the light will come on in the chamber.

Molex Pin #	Function	Wire Color
1	Not Used	
2	Operate	White
3	+28 Volts	Red

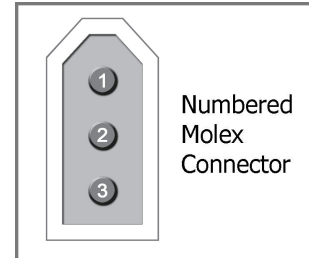

APPENDIX A | CONTACT INFORMATION

Please contact MED Associates, Inc. for information regarding any of our products.

Visit our website at www.med-associates.com for contact information.

For technical questions, email support@med-associates.com.